




Double Safe Waste Chute SOP

1. Carefully remove inner waste chute bag, wrap and twist, then seal with rubber band and discard into the waste chute.




Carefully Remove Inner Waste Chute Bag, Wrap and Twist, Then Seal With Rubber Band, and Discard into Waste Chute.

2. Take new bag and wrap rubber band around open end. Place bag, facing into enclosure, on the mount.

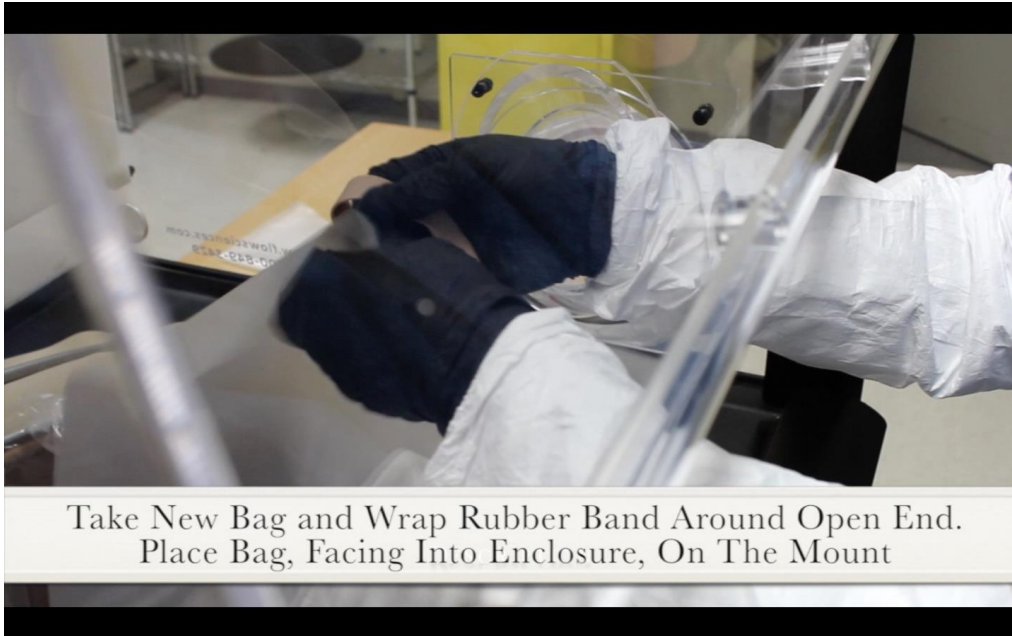


3. Carefully remove sleeve protectors, then outer glove layer, and place inside of the enclosure.

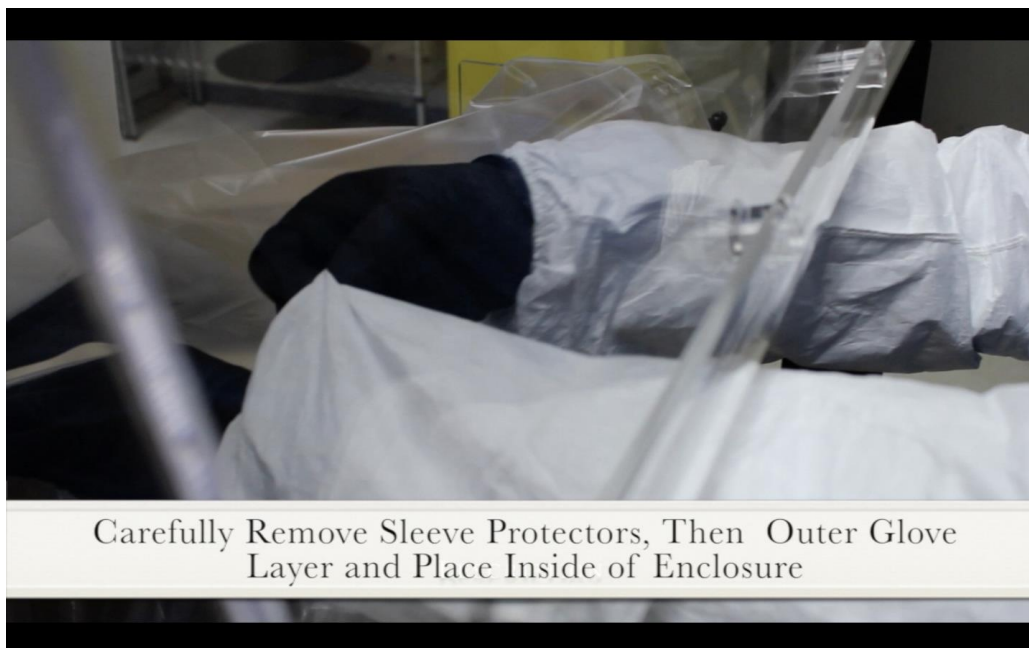

4. Remove outer bag, wrap and twist, and seal with rubber band. Then discard into trash can.



5. Reach through the new bag, and grab gloves and sleeve protectors from inside. Pull all through to outside of enclosure.

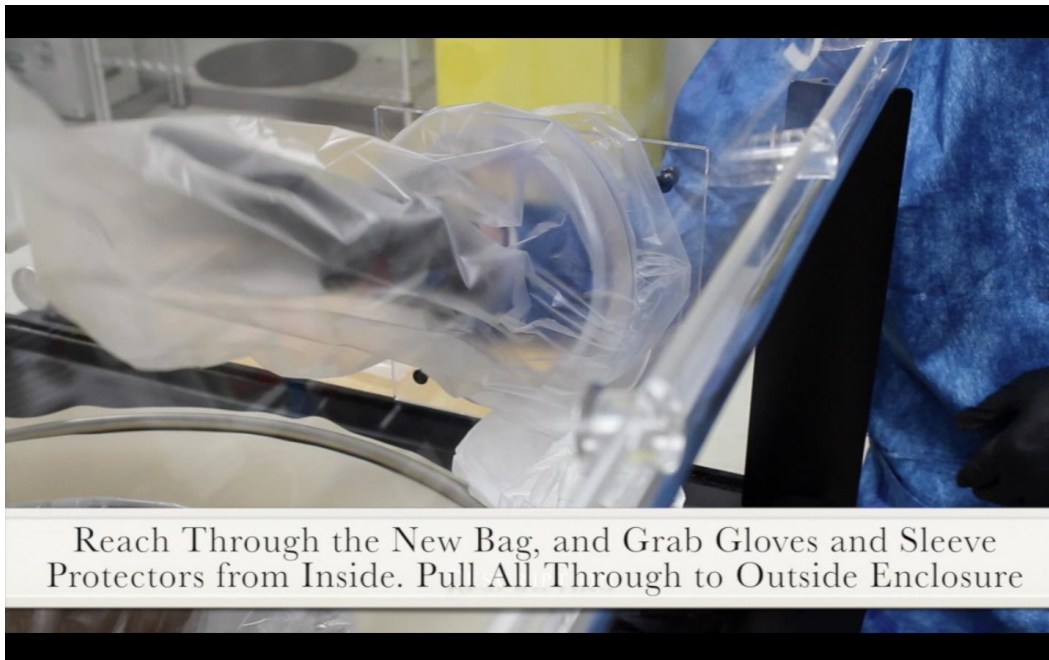

6. Place new bag onto waste chute from outside. Place OVER the bag that was pulled through.



7. Remove glove layer and discard into trash.



You have now completed the Double Safe Waste Chute
change out process.